

VisionGraft®

THE CLEAR CHOICE®

A sterile, non-immunogenic, cross-linked cornea, successfully used in more than 70,000 ocular surgeries.

Applications for Cornea

- Treats corneal melts, ulcers, and perforations
- Anterior Lamellar Keratoplasty
- Tectonic/therapeutic PKP
- Keratoprosthesis ring



Engineered Precisely to Meet Your Clinical Demands

①

Cryogenically treated
to inactivate cells

②

Processed specifically
for a range of
corneal surgeries

③

Terminally sterilized
with precision gamma
irradiation to a
sterility assurance
level of 10^{-6}

④

Stored in a
proprietary medium
that is largely
albumin based

For more information or to order,
please contact us at 800-858-2020
or visit www.CorneaGen.com

CorneaGen™

Tissue Codes

Code	Description	Shape
C0101AL	Whole Cornea with Scleral Rim	
C0100AL	Whole Cornea without Scleral Rim	
C0300AL	Split Thickness Full Moon	
C0400AL	K-Pro Ring	

VisionGraft is sterile, ready to use, and remains as clear as non-irradiated corneas.¹

NON-IMMUNOGENIC

Reduces risk of infection for patients.

- No adverse reactions reported
- All cells removed during processing
- Eliminates risk of bacterial or fungal disease

READILY AVAILABLE

Stable at room temperature for up to 2 years.

- Easy to stock for emergencies
- Variety of sizes for various surgical applications
- Engineered from medically eligible corneas that are not ideal for traditional transplant

CUSTOM SIZE OPTIONS


CorneaGen can process VisionGrafts to your specifications.

- Special orders can be completed quickly
- Collaborate to explore new applications for use


DURABLE AND DEPENDABLE

Biomechanically comparable to fresh corneas.

- Histology shows structure is intact post-irradiation²
- Suture strength comparable to fresh corneas²
- Tissue quality comparable to fresh corneas

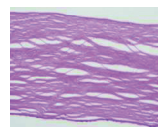

Irradiated Cornea

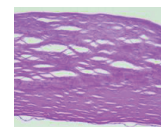

Fresh Cornea

REFERENCES

1. Chae JJ, Choi JS, Stark WJ, Elisseff J. Extracellular matrix characterization of the acellular gamma-irradiated cornea. Poster presented at: The Association for Research in Vision and Ophthalmology Annual Meeting May 5-9, 2013; Seattle, WA. DOI C0176.
2. Daoud YJ; Smith R; Smith T; Akpek EK; Ward DE; Stark WJ. The intraoperative outcomes of gamma-irradiated corneas in corneal and glaucoma patch surgery. *Cornea* 2011; 30(12): 1387-1391.

Exclusively processed and distributed by:

CorneaGen™

101 N. Chestnut Street
Ste 303
Winston-Salem, NC 27101 USA

800-858-2020
Bal.Tissue@CorneaGen.com
Products@CorneaGen.com